

CAPTURE BT FUMING CHAMBER ACCESSORIES



1 ALUMINUM FUMING DISHES

These dishes are ideal for Cyanoacrylate fuming. A small amount of cyanoacrylate placed inside the dish and placed on a hot plate or in a chamber will cause rapid fuming of the cyanoacrylate.

A-2558D 2.75 x 5/8" 2.7 oz capacity 100/pk

3 RHODAMINE 6G DYE STAIN

Rhodamine 6G is one of the best dyes to use with cyanoacrylate developed latents and a laser or forensic light source. Available as a premix liquid for immediate use, or as a premeasured concentrate to be mixed with methanol or distilled water.

A-26481 25 gm
A-PF045 1 Liter Premix

2 CYANOACRYLATE

Liquid cyanoacrylate designed specifically for use in fingerprint development. Liquid glue is placed in a small metal container and placed on a hot plate with a source of humidity. The cyanoacrylate (CA) will vaporize into cyanoacrylate esters that will adhere to latent prints, making them a light, chalky white. The resultant latent prints can then be photographed or further processed with powder, dyes or stains.



A-3202B 2 oz
A-2616 16 oz

4 CYANO-FORCE™ FUMING CHAMBER WIPES

CYANO-FORCE™ Citrus Fuming chamber wipes combines the strength of a natural citrus solvent with a highly absorbent, abrasive, nonscratching towel for a one-step cyanoacrylate removing solution. It emulsifies and traps the glue, eliminating a chance of redepositing on the cleaned surface. Provides controlled application and a safer alternative to using harsh solvents.

A-90630 30 Wipes, 10.5" x 12.25"

CAPTURE BT FUMING CHAMBER ACCESSORIES

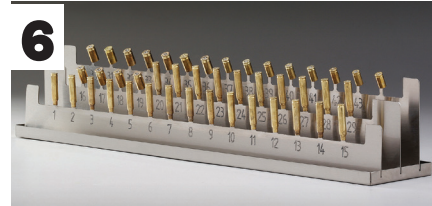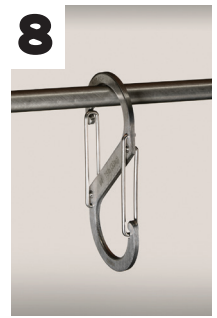

1. Capture BT Mobile Stand. Powder-coated steel frame with 5" diameter, toelocking, polyurethane casters. Manually adjustable height from 21" to 28" h (53 to 71 cm). Shelves and compartment hold supplies and accessories. Overall width x depth: 37.6" x 34.5" (96 x 88 cm). Shipping weight 125.0 lbs. (57 kg)
A-3182900 Mobile Cart

2. Wire Shelf Kits. Each shelf fits onto CAPture BT sidewall brackets. The CAPture BT can hold up to four wire shelves. Shelf weight limit: 25 lbs. (11 kg). Stainless steel. Shelf dimensions: 26" w x 23" d (66 x 58 cm).
A-3181511LC Single Wire Shelf

3. Wire Half Shelf Kit. Includes two wire half-depth shelves. Each shelf fits onto CAPture BT sidewall brackets. The CAPture BT can hold up to eight half-depth shelves. Stainless steel. Shelf weight limit: 15 lbs. (7 kg). Shelf dimensions: 26" w x 12" d (66 x 30 cm). Shipping weight 5 lbs. (2 kg)
A-3181522LC Half Wire Shelf Kit - 2 Shelves

4. Hanging Rod Kit. Includes four stainless steel rods, 26.9" long x 3/8" diameter (68 x 1 cm). Rod weight limit: 15 lbs. (7 kg). Four each rods are included with each CAPture BT. Shipping weight 5.0 lbs. (2 kg)
A-3184000LC Hanging Rod Kit

5. Perforated Shelf Kits. Each shelf fits onto CAPture BT sidewall brackets. The CAPture BT can hold up to three perforated shelves. Stainless steel. Shelf dimensions: 26" w x 23" d (66 x 58 cm). Shelf weight limit: 25 lbs. (11 kg).
A-3181611LC Perforated Shelf Kit - 1 Shelf

6. Casing Holder Kit. Provides 45 numbered standoffs for upside down casing placement. For casings .32 caliber or larger. Shipping weight 10 lbs. (5 kg)
A-3181700LC Casing Holder Kit

7. Long Gun Holder Kit. Provides five U-shaped cradles for long guns to be held securely. Holder locks into existing side wall brackets. CAPture BT can accommodate up to 15 long guns (three Holder Kits). Stainless steel. Shipping weight 4 lbs. (2 kg)
A-3184100LC Long Gun Holder Kit

8. Large Clip Kit. Includes ten large clips. Each clip is 3.5" x 1.5" x approximately 0.3" thick (9 x 4 x 0.8 cm). For easy hanging of larger items, such as hand guns. Stainless steel. Shipping weight 2 lbs. (1 kg)
A-3184200LC Large Clip Kit

9. Mini Clip Kit. Includes ten small, alligator-style clips. Each clip is 1.2" x 0.3" (3 x 0.8 cm). For easy hanging of smaller, lightweight items, such as plastic bags. Shipping weight 1 lb. (0.5 kg)
A-3184300LC Mini Clip Kit

FILTERS

10. Carbon, Main, Pre-filter (10/pk). Replace every 50 cycles.
A-3181400LC 10/pk

11. Replacement Humidifier Filter. Replace every 500 cycles.
A-3181401LC

12. Replacement Carbon Filter, Exhaust. Replace every 500 cycles.
A-3186500LC

13. Replacement Complete Filter Kit Provides all filters necessary for approx. 500 cycles.
A-3185900LC